





When you consider comfort, quality and performance for your healthcare mattress, Simmons™ Clinical Care is at the top of the list. Simmons™, maker of consumer mattress icons such as BeautyRest® and ComforPedic®, introduces a full line of therapeutic sleep surfaces that balance comfort with clinical outcomes. Providing a good night's rest is as important in a healthcare setting as it is in your home. Simmons™ Clinical Care mattresses manage pressure, shear and heat by using the latest design enhancements and exclusive materials to ensure resident comfort. A full line of static foam core products with varying levels of immersion help with pressure redistribution. By modifying the surface we reduce shear and allow for heat to be pulled away from the body. The Simmons™ Clinical Care Plus, Sensor Plus and Max provide the latest in active support surface therapy. Alternating Pressure, Low Air Loss and True Low Air Loss therapies provide the next level of pressure redistribution and heat/moisture control for your residents at highest risk.



S600 Series Foam Mattress



With Firm Side Bolsters



With 3" Raised Firm Side Bolsters



Cutaway Closeup


- Protection for higher risk patients with pressure ulcers and high-risk heels
- Provides pressure reduction and redistribution
- Firm side perimeter construction with optional 3" raised perimeter bolsters (RB) featuring ingress-egress cutout
- Tri-laminate construction with heel slope
- Fluid-resistant, low friction, anti-microbial top cover with vinyl bottom cover and full zipper
- Designed to meet the requirements of HCPCS Code E0184
- Maximum weight capacity, EVENLY DISTRIBUTED:
35" wide models: 450 lb
42" wide models: 500 lb
- Warranty: mattress warranty five-year limited; cover warranty one-year limited

Meets Entrapment Guidelines when used with Zenith9000APS, Zenith9000, Zenith7000, Zenith5000, Matrix II, Liberty, and B700



| Dimensions | With Firm Side Bolsters | With 3" Raised Firm Side Bolsters | Qty |
|----------------|-------------------------|-----------------------------------|------|
| 35" x 76" x 6" | S600B-3576 | S600RB-3576 | 1 ea |
| 35" x 80" x 6" | S600B-3580 | S600RB-3580 | 1 ea |
| 35" x 84" x 6" | S600B-3584 | S600RB-3584 | 1 ea |
| 42" x 76" x 7" | S600B-4276 | S600RB-4276 | 1 ea |
| 42" x 80" x 7" | S600B-4280 | S600RB-4280 | 1 ea |
| 42" x 84" x 7" | S600B-4284 | S600RB-4284 | 1 ea |

S500 Series Foam Mattress

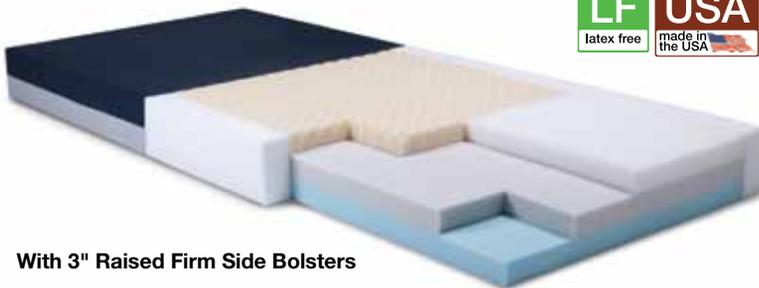

With 3" Raised Firm Side Bolsters



With Firm Side Bolsters



Cutaway Closeup

- Preferred combination of pressure redistribution and support
- Modified convoluted topper enhances immersion characteristics and increases airflow
- Tri-laminate construction with soft heel-zone
- Firm bolster with optional 3" raised bolster
- Fluid-resistant, low friction, anti-microbial top cover with vinyl bottom cover and full zipper
- Designed to meet the requirements of HCPCS Code E0184
- Maximum weight capacity, EVENLY DISTRIBUTED:
35" wide models: 400 lb
42" wide models: 450 lb
- Meets 16 CFR 1633
- Warranty: mattress warranty five-year limited; cover warranty one-year limited

Meets Entrapment Guidelines when used with Zenith9000APS, Zenith9000, Zenith7000, Zenith5000, Matrix II, Liberty, and B700



| Dimensions | With Firm Side Bolsters | With 3" Raised Firm Side Bolsters | Qty |
|----------------|-------------------------|-----------------------------------|------|
| 35" x 76" x 6" | S500B-3576 | S500RB-3576 | 1 ea |
| 35" x 80" x 6" | S500B-3580 | S500RB-3580 | 1 ea |
| 35" x 84" x 6" | S500B-3584 | S500RB-3584 | 1 ea |
| 42" x 76" x 7" | S500B-4276 | S500RB-4276 | 1 ea |
| 42" x 80" x 7" | S500B-4280 | S500RB-4280 | 1 ea |
| 42" x 84" x 7" | S500B-4284 | S500RB-4284 | 1 ea |

S400 Series Foam Mattress



With Firm Side Bolsters



With 3" Raised Firm Side Bolsters



Cutaway Closeup

- Medium acuity mattress for patients with existing low grade ulcers
- Provides pressure reduction and redistribution
- Firm side perimeter construction with optional 3" raised perimeter bolsters (RB) featuring ingress-egress cutout
- Three-Layer construction with special heel section
- Convoluted *WAVE* foam top surface enhances air flow to reduce moisture buildup and skin maceration
- Fluid-resistant, low friction, anti-microbial top cover with vinyl bottom cover and full zipper
- Meets 16 CFR 1633 requirements
- Designed to meet the requirements of HCPCS Code E0184
- 350 lb maximum weight capacity, **EVENLY DISTRIBUTED**
- Warranty: mattress warranty four-year limited; cover warranty one-year limited

| Dimensions | With Firm Side Bolsters | With 3" Raised Firm Side Bolsters | Qty |
|----------------|-------------------------|-----------------------------------|------|
| 35" x 76" x 6" | S400B-3576 | S400RB-3576 | 1 ea |
| 35" x 80" x 6" | S400B-3580 | S400RB-3580 | 1 ea |
| 35" x 84" x 6" | S400B-3584 | S400RB-3584 | 1 ea |
| 42" x 76" x 6" | S400B-4276 | S400RB-4276 | 1 ea |
| 42" x 80" x 6" | S400B-4280 | S400RB-4280 | 1 ea |
| 42" x 84" x 6" | S400B-4284 | S400RB-4284 | 1 ea |



Meets Entrapment Guidelines when used with Zenith9000APS, Zenith9000, Zenith7000, Zenith5000, Matrix II, Liberty, and B700

S300 Series Foam Mattress



With Firm Side Bolsters



With 3" Raised Firm Side Bolsters



Cutaway Closeup

- Low acuity mattress — first line of defense against pressure ulcers for “at-risk” patients
- Provides pressure reduction and redistribution
- Firm side perimeter construction with optional 3" raised perimeter bolsters (RB) featuring ingress-egress cutout
- Three-Layer construction with comfort heel section
- “Flat-Foam” top surface provides comfort and pressure reduction
- Fluid-resistant, low friction, anti-microbial top cover with vinyl bottom cover and full zipper
- Meets 16 CFR 1633 requirements
- Designed to meet the requirements of HCPCS Code E0184
- 300 lb maximum weight capacity, **EVENLY DISTRIBUTED**
- Warranty: mattress warranty three-year limited; cover warranty one-year limited

| Dimensions | With Firm Side Bolsters | With 3" Raised Firm Side Bolsters | Qty |
|----------------|-------------------------|-----------------------------------|------|
| 35" x 76" x 6" | S300B-3576 | S300RB-3576 | 1 ea |
| 35" x 80" x 6" | S300B-3580 | S300RB-3580 | 1 ea |
| 35" x 84" x 6" | S300B-3584 | S300RB-3584 | 1 ea |
| 42" x 76" x 6" | S300B-4276 | S300RB-4276 | 1 ea |
| 42" x 80" x 6" | S300B-4280 | S300RB-4280 | 1 ea |
| 42" x 84" x 6" | S300B-4284 | S300RB-4284 | 1 ea |



Meets Entrapment Guidelines when used with Zenith9000APS, Zenith9000, Zenith7000, Zenith5000, Matrix II, Liberty, and B700

Clinical Max Low Air Loss with Alternating Pressure System



The Clinical Max is a low air loss mattress replacement system that combines TRUE low air loss and dynamic alternating pressure — you no longer have to sacrifice one therapy for another. With its high air flow and water vapor permeable cover, the system can aid in the prevention of maceration and moisture build up. With various therapy features including adjustable cycle time, pulsation therapy, static mode and seat inflation, the Clinical Max provides optimal pressure relief.

- TRUE Low Air Loss technology wicks moisture away from patient
- Air output: > 1340 LPM
- 10 comfort levels for optimizing therapeutic support and patient comfort
- 2" pocket for optional foam base
- Adjustable alternating cycle time meets the treatment needs for each individual
- Oscillated pulsation mode helps improve circulation
- Seat inflation mode prevents patient from bottoming out while seated
- Safety alarms with audio and visual indicators to alert caregivers
- Smart vacuum for rapid deflation and quick storage
- MaxFirm for easy patient transfers and nursing procedures
- Designed to meet the requirements of HCPCS Code E0277
- 440 lb maximum weight capacity, EVENLY DISTRIBUTED
- Two-year limited warranty on control unit; one-year limited warranty on mattress

770000S 79" x 36" x 8"
(10" with optional foam base) 1 ea

763400F 1 ea Optional Foam Base

Clinical Sensor Plus Alternating Pressure Mattress System with Active Sensor Technology



- Intelligent, alternating pressure relief support surface that uses proprietary Auto-Adjusting technology, used to automatically configure air pressure based on patient weight, relieving caregivers of any manual setup and assuring maximum comfort for the patient
- Group II Mattress System suitable for use in the prevention and treatment of pressure ulcers
- Touch pad adjustment with LED indicators
- Cycle Time: 10 minutes
- Fine tuning function to manually adjust the air pressure settings
- Four therapy settings: Alternating Pressure, Static, Alternating + Seat Inflation and Static + Seat Inflation
- Lockable control panel to prevent unauthorized adjustment changes
- Cell-on-cell design ensures that the mattress will hold air in the lower cells in the event of a power failure, preventing the user from bottoming out
- Low pressure and power failure LED indicators and audible alarms
- Quilted nylon polyurethane cover is vapor permeable and fluid resistant
- Quick CPR deflation
- 400 lb maximum weight capacity, EVENLY DISTRIBUTED
- Designed to meet the requirements of HCPCS Code E0277
- One-year limited warranty

765000S 78.7" x 35.4" x 8" 1 ea

Clinical Plus Alternating Pressure / Low Air Loss Mattress



This ideal, portable air flotation system, intended to help and reduce the incidence of pressure ulcers while optimizing patient comfort, offers both alternating pressure and low air loss. Quick CPR deflation is accomplished by simply disconnecting quick-connect hose couplings from control unit.

- Suitable for use in the prevention and treatment of stages 1-3 pressure ulcers
- Touch pad adjustment with LED indicators
- Ten comfort settings for maximum patient comfort and adjustability
- Four therapy settings including autofirm for inflation, alternating pressure, static, or seat inflation
- Cell-on-cell design ensures that mattress will hold air in the lower cells in the event of a power failure, prevent the user from bottoming out
- Low pressure and power failure LED indicators and audible alarms
- Pocket for optional 2" foam base
- Quilted nylon/polyurethane cover
- Air input > 8 LPM
- Cycle time 10, 15, 20, or 25 min
- Electrical rating: 120 VAC/60Hz
- UL60601-1 and CSA C22-2 listed
- Designed to meet the requirements of HCPCS Code E0277
- 400 lb maximum weight capacity, EVENLY DISTRIBUTED
- One-year limited warranty

760000S 1 ea

763400F 1 ea Optional Foam Base



BASIC AMERICAN
MEDICAL PRODUCTS

View these and more products at www.grahamfield.com

770-368-4700


GF Health Products, Inc. is not responsible for typographical errors. All illustrations, specifications, packaging and warranties contained in this literature are based on the latest product information available at the time of printing. The most current product information can be found online at www.grahamfield.com © 2013, GF Health Products, Inc.

SIMMONS Clinical Care products manufactured by GF Health Products, Inc. under license from the trademark owner, Dreamwell, Ltd., a subsidiary of Simmons Bedding Company, Atlanta, Georgia.

SIMMONS is a trademark of Dreamwell, Ltd. a subsidiary of Simmons Bedding Company. © 2013 Dreamwell, Ltd. All rights reserved.

GF1200146RevA13

Lightest Gel Overlay

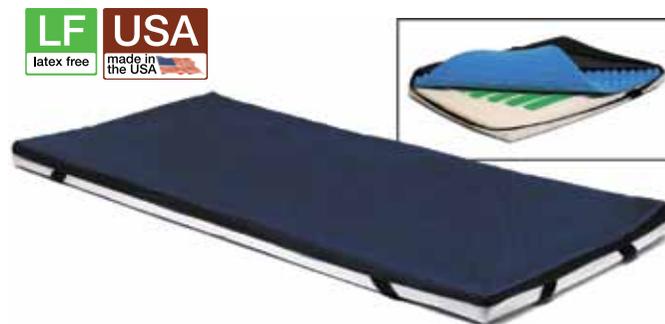

- A new standard in pressure management
- Extremely lightweight for easy handling and lower shipping cost
- Mattress overlay incorporates a solid gel patterned layer that is inherently visco-elastic, and works in responding to pressure by deforming slightly to disperse critical weight over a large surface area
- Nonmigrating gel is engineered into a high-density 2" solid foam base with specially-designed 1" convoluted top to provide maximum comfort and pressure management
- Overlay reduces body interface pressure and provides comfort at critical points over a wide range of body shapes and sizes
- Nonaqueous gel will not migrate or pool, making the overlay effective in any bed position
- Overlay delivers same consistent results over lifetime of product
- Fluid-resistant, antibacterial, antifungal cover meets CAL 117
- Two-year limited warranty

• 275 lb maximum weight capacity, EVENLY DISTRIBUTED
103680S 80" x 36" x 3" 1 ea

• 275 lb maximum weight capacity, EVENLY DISTRIBUTED
103684S 84" x 36" x 3" 1 ea

• 400 lb maximum weight capacity, EVENLY DISTRIBUTED
104280S 80" x 42" x 4" 1 ea