

Zenith 7000APS

Advanced Positioning System

 **Basic
American**™



Zenith 7000APS

Advanced Positioning System

The advanced positioning of the **Zenith7000APS** allows the bed to go into numerous positions for:

- Positively influencing clinically driven outcomes through multiple positioning options
- Accommodating post-surgical rehab positioning to facilitate resident comfort
- Reduces caregiver strain and potential injuries to caregiver and residents caused by manually repositioning residents



Zenith7200APS

SilverSolutions within the powder-coating provides antimicrobial surface protection.

SilverSolutions inhibits the growth of organisms that can detract from the appearance and coating functionality:


- Inhibits the growth of bacterial odors
- Resists stains caused by mildew
- Resists stains caused by mold

SilverSolutions preserves the cured coating during the service life of the bed.

GF Health Products, Inc. makes no independent representations or empirical statements about SilverSolutions.

This information is reprinted from: TCI Powder Coatings, EPA Guidance for EPA Registration 40 CFR 152.25



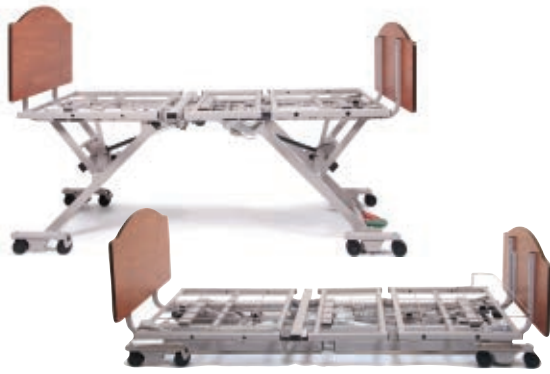

Designed for comfort and efficient operation

Height travel range: 8.75" to 27"

Straight vertical rise protects walls and maintains footprint.
Timesaver Technology™ - fast motors save time.

Strong, no-sag grid deck sleep surface with sturdy cross braces.

Safe work load: 500 pounds*



Sealed ball bearings in over 20 locations ensure maintenance-free, continuous, silent operation.

Bed provides 360° mobility at ANY height, ideal for cleaning under the bed.

Head and foot boards are securely mounted with tall, tubular mounting reduces breakage, creating a built-in push bar.

Bed locks and rolls at any height, reducing caregiver back strain.

Tool-less assembly reduces set-up time.



Multiple positioning such as lounge chair and recliner

style functionality provides comfort and relaxation to the resident and also allows for positioning of resident closer to the caregiver for proper care.

The **Cardiac Chair** is primarily designed for residents with heart and respiratory problems to whom a semi-erect position brings great relief. The chair can be made to serve as a fully reclined bed when needed, but can also be adjusted to form a chair.

This "cardiac position" helps to elevate the patient from the waist up, which provides relief to the lungs, increases circulation, and assists the patient in coming from a fully prone position to a sitting position without causing harm or undue strain.

Standard features of the Zenith7200APS



Maximum underbed clearance for use with overbed tables and lifts.



Pendant Hand Control and Holster

Easy operation contoured pendant is shock-resistant, splash-proof and durable. On the Zenith7200APS, pendant connects to either side of the bed to adapt to varying room layouts. Holster conveniently secures to bedrails and assist devices for easy access.



Intuitive two-pedal locking system

is centrally positioned and easily visible from a distance to ensure that the bed is locked. Pedal locks are readily accessible at ALL bed heights and widths, even in the lowest position. Bed rolls at any height reducing caregiver back strain.



• Adjustable wall saver

positions (76", 80", and for use with a trapeze) provide the smallest footprint at each length, while protecting the wall (included with bed).




♦ Trendelenburg and Reverse Trendelenburg

positioning, ideal for rehab and subacute, provides therapeutic benefits that can aid in circulation and ease body discomfort. Trendelenburg also promotes gravity assist resident self positioning.

♦ Bed automatically pauses when it returns to level from Trendelenburg, Reverse Trendelenburg or chair position.

♦ Auto CPR functionality allows for ideal CPR positioning in emergency situations.

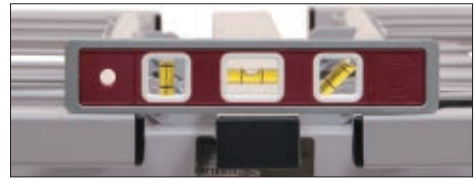

Synchronized high-low motors

keep the bed level for proper therapeutic positioning. When the head position is lower than the heart, lung volume may be compromised, causing airway problems, making it more difficult to breathe. Uneven sleep surfaces may cause a risk of brain swelling and reduced systolic and diastolic pressures.



Grounded bed with right-angle hospital-grade plug is a standard feature.

• Adjustable mattress retainer

securely supports mattress position for 76" or 80" length (included with bed).



Safety inspired **embedded staff control** in footboard

- Angled for better visibility.
- Two-stage lockouts – lockout of pendant use and lockout of embedded staff control.
- Quick lockout buttons lock multiple functions simultaneously.
- Industry Exclusive **SafeMode** can automatically lock out the high/low feature, if necessary.

Directional mechanism on front casters is accessible at all heights for ease in rolling the bed.




Optional Half Length Assist Device and Assist Bars provide a sturdy and secure hand hold to assist residents in and out of bed. Assist devices are removable without tools and meet FDA Entrapment Guidelines.




Warm-to-the-touch Half Length Assist Device
ZA84200

Warm-to-the-touch half-length assist device and pivot assist bars:

- Move out of the way for lateral transfers.
- Pivot away from egress area.
- Provide toe clearance even in lowest position to prevent staff injury.
- Hand pendant mounts facing inward or outward and away from hand grab area.

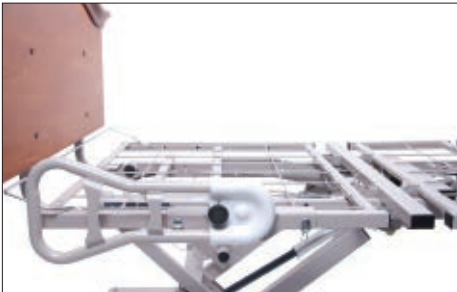

Pivot Assist Bar ZA85400 (tool-less)
ZA87900 (bolt-on)

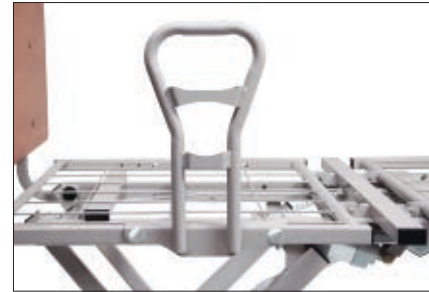

Fixed Assist Bar ZA85000 (tool-less)
ZA87800 (bolt-on)



Extension Kit


Transform your bed to 84" with this optional 4-inch sleep deck extension.

ZA85350 Standard 35" width
ZA85300 39" or 42" widths

Battery Backup

Provides the ability to perform bed functions in the event of power outage.

999-0711-003 Battery backup
999-0711-002 Charger

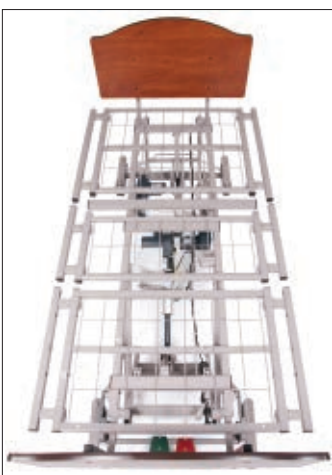

Wide Deck Kit

Optional wide kit adjusts the bed sleep deck to 39" or 42" width for greater resident comfort. (Head and foot boards and staff control electronics not included.)

ZA858000 39" wide deck kit
ZA851000 42" wide deck kit



May be attached with bolts or quick release pins.

Trapeze Adapter

Optional trapeze adapter allows for easy attachment of trapeze system to headboard within easy reach for residents for repositioning themselves.



ZA79000 Trapeze adapter
ZA78100 Trapeze bar, chain, and triangle

APS78174 Zenith7200APS with pedal lock

U.S. Patents

Height Travel Range.....	8.75" to 27"
Safe Work Load*.....	500 lbs.
Mattress Deck Length.....	76"/80"
Mattress Deck Width.....	35"
Overall Bed Length (76"/80")	81.75"
Overall Bed Width.....	39" with assist bar
Overall Bed Width.....	40.75" with half length assist devices
Maximum back angle.....	70° (Chair position 84°)
Maximum knee angle.....	30°



- 15 years on frame
- **5 years on motors and control box**
- 3 years on other electronics

For more information on the Zenith Series, please visit us online at www.grahamfield.com or call Basic American today.

Highest Safety and Quality Standards Unrestricted UL Listing

- **Bed approved as a transport device because beds pass the threshold clearance/impact testing at ALL heights by 20mm (roughly 3/4")**
- No pinch points within 200mm of bed perimeter
- Splash-proof electronics
- Load capacity and tip-stability tested and approved



Information contained herein is subject to change.



GF1500004 RevB16