



EXTENDED CARE BEDS



We proudly make all of our beds in the USA to the highest safety standards:⁴

- FDA Recognized Standard: ANSI/AAMI STD ES60601-1
- Health Canada Recognized Standard: CAN/CSA C22.2 No. 60601-1 (IEC 60601-1:2012-Edition 3.1)
- Includes International Standards: IEC 60601-1, IEC 60601-1-2, IEC 60601-2-52

All Basic American Extended Care Beds Are

- Certified transport device - passes threshold clearance at ALL heights - ideal for evacuation
- Silver Solutions⁸ powder-coating provides antimicrobial surface protection and deters the growth of odor-causing bacteria
- No holes on top surface of tubes that could harbor fluids/bacteria
- Grounded bed with right-angle hospital grade plug for safety
- Adjustable Wallsaver protects the wall with smallest footprint
- Patented Hi-Lo mechanism adds life to the motors
- Sturdy, no-sag, grid sleep surface with cross braces
- Bed provides 360 mobility for cleaning, etc.
- Proudly Made in the USA
- Assist devices available*
- Tool-less assembly



Adjustable Wallsaver

- Maximum wall protection
- Adjustable clearance lengths
- Standard on all beds

Optional
Zenith 7200
Pendant Shown



Slide-W-I-D-E[®]
Models Only

Mattress Retainers⁹

- Adjustable for multiple mattress lengths
- Standard on all beds



All Other
Models

Zenith & Matrix Series Bed Standard Features

- Advanced Positioning System (APS) includes One-Button Comfort/Cardiac Chair Position (optional on Matrix 7200).
- Straight-Line Hi-Lo Travel with Patent-Awarded Design Technology
- 'Plug & Play' Digital Electronics enable future feature additions
- Smart Technology
- Industry-leading, energy-saving electronics consume up to 66% less power than transformer electronics when idle
- Sealed ball bearings joints provide continuous, maintenance-free, silent operation
- Intuitive two-pedal locking system - centrally located, easily visible, and accessible at all heights and widths
- Tool-less assembly

Zenith & Matrix Series Bed Configurations

- Embedded staff control with safety-inspired, two-stage lockouts enabling caregivers to lock-out features to the resident
- Standard or Advanced Positioning System (APS)
- APS Staff Controls are optional. (APS Staff Control includes One-Button Comfort/Cardiac Chair and Trendelenburg/Reverse Trendelenburg Positions)
- Hand Pendant to control Trendelenburg/Reverse Trendelenburg (sold separately)
- Accessories and some options are sold separately*



Intuitive Two-Pedal Locking System

- Centrally positioned and easily visible from a distance ensures the caregiver that bed is locked at a glance.
- Bed rolls at any height increasing caregiver safety
- Readily accessible at ALL bed widths and heights (even lowest)
- Select Option at time of purchase



Zenith 9200 SERIES

9200 WITH APS SHOWN WITH
Slide-W-I-D-E® deck option

Zenith

Zenith Series Standard Features

- Straight-Line Hi-Lo Travel with Patent-Awarded Design Technology
- Mounted Head and Footboards create a built-in push bar and reduce breakage




Zenith 7200 SERIES

7200 & 7200 WITH APS

OPTIONAL WITH THE ZENITH SERIES

Slide-W-I-D-E®

EXPANDABLE/RETRACTABLE SLEEP DECK

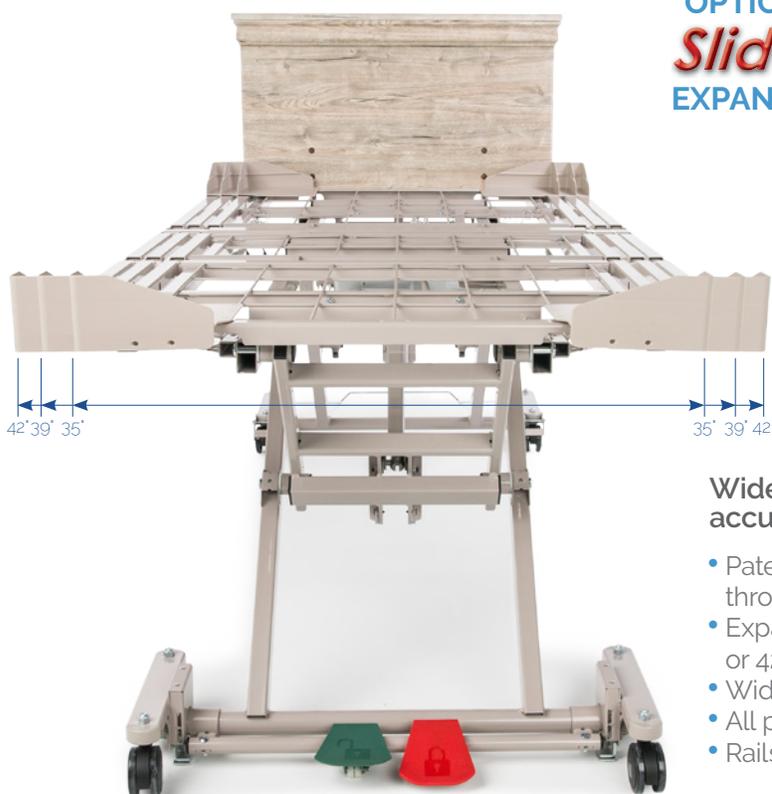

Quick-release
ball-detent pins
attached to beds
with lanyards



Wider sleep surfaces may help reduce falls by residents accustomed to sleeping on common consumer products.⁷

- Patented design allows for quick adjustment, ideal for moving through doorways to assist in emergency evacuations
- Expandable/retractable deck quickly and easily adjusts to 35, 39, or 42 inches with two quick-release ball-detent pins on each side.
- Wide head and footboards available in 42 inch width
- All pins are attached to the bed with lanyards, eliminating lost parts
- Rails and assist devices are quick-release without tools



Matrix 6200 APS SHOWN
with intuitive pedal-lock option
Raised Head Position

Matrix

Matrix Series Standard Features

- Advanced Positioning System (APS) includes One-Button Comfort/Cardiac Chair Position (optional on Zenith 7200).
- Straight-Line Hi-Lo Travel with Patent-Awarded Design Technology

Matrix 6200 SERIES 6200 WITH APS

Advanced Positioning System (APS)

Advanced positioning allows the bed to go into numerous positions

- Positively influencing clinically driven outcomes through multiple positioning options
- Accommodating post-surgical rehab positioning to facilitate resident comfort
- May reduce caregiver strain and potential injuries
- Built-In Hi-Lo Lockout Feature (on APS Hand Controls)



Reverse Trendelenburg



Trendelenburg Position



Zenith 9200 APS

SHOWN
with Advanced Positioning System, Slide-W-I-D-E®
deck option, and intuitive pedal-lock option
Comfort/Cardiac Chair Position



Liberty SHOWN
with intuitive pedal-lock option
Supine Position



Liberty Bed Standard Features

- Desired features at affordable price point
- Four Locking casters

Liberty Bed Available with

- Optional one-step red/green pedal lock to maximize your security (*as shown*)
- Optional 4-inch Total-lock casters increase under-bed clearance

Overbed Table Options





Rectangle



Raised Panel



Mission




Suwanee³



Center Crown




Crescent with Cap³

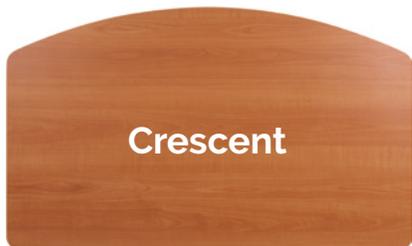

Crescent

3DL material and custom finishes available for all board styles upon request - costs may vary

Counter-Rotating Assist Device

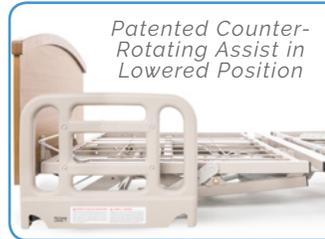


Patented Counter-Rotating Assist in Lowered Position



Raised Position

Provides a sturdy and secure hand hold to assist residents in and out of bed yet easily swivels out of the way when desired.

- Meets FDA Entrapment Guidelines
- Patented Counter-Rotating Assist Devices
- Removable without tools
- Moves out of the way for lateral transfers
- Pivots away from egress area
- Provides toe clearance in lowest position to help prevent staff injury
- Hand pendant mounts facing inward or outward and away from hand grab area

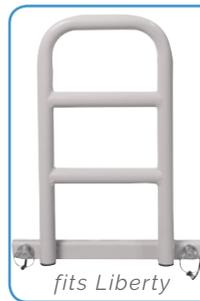

fits Liberty



fits Matrix & Zenith

Pivot Assist Bars



Pivot Assist Device - Rotated Down, Foot



Pivot Assist Device - Up-Right Position



Pivot Assist Device - Rotated Down, Head

Fixed Assist Devices

Provides a sturdy and secure hand hold to assist residents in and out of bed.

- Moves out of the way for lateral transfers
- Tool-free assembly and removal



Bed Cart

- Ideal for ease of storage.
- Allows beds to fit through narrow doorways.
- Designed for One-Person Mounting and maneuverability
- Compatible with all Basic American models



USB Power Supply

- Charges most devices
- Compatible with all Basic American models



Side Wall Bumper

- Maximizes wall protection
- Compatible with all Zenith and Matrix models

Advanced Positioning Pendant

- Upgrade from Standard to APS on Zenith 7200
- One-Button Comfort/Cardiac Chair Position
- Hi-Lo Lockout Feature





Board Mounted IV Socket

- Mounts with 2 screws
- Easy to install
- IV pole sold separately
- Compatible with all Basic American models²



Surface Length Extension Kit

- Extends standard length bed frame from 84" to 88". Complete with solid pan extension, mattress guard, and installation hardware.
- Compatible with all Basic American models²
- Easy to install

UnderBed Light

- Helps reduce fall risks
- Eliminates bright light intrusion
- Multi-level light intensity
- Compatible with all Basic American models



Embedded Staff Control

- Standard Functionality and Lockouts, or Advanced Positioning (APS).
- Angled-Up Display for Improved Visibility
- Two-Stage lockouts deter unwanted changes
- Mix or Match with Plug & Play design
- Upgrade anytime²




Plastic Pan Deck Kit

- Compatible with all Zenith and Matrix models²



Edema Foot Ratchet Kit

- Elevated foot rest.
- Standard on Zenith 9200
- Optional on Matrix 6200 and Zenith 7200

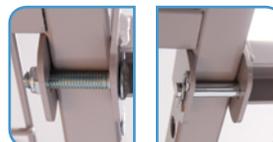

Bolt & Quick Release Pin shown on Wide Kit



Pin-on Wide Kit

- Compatible with all 35" wide Zenith and Matrix models
- Available in 39" or 42" width kits
- Anytime add-on²



Trapeze and Adapter

- Pivots out of the way when not in use
- Compatible with all Zenith and Matrix models
- Rigid Design

Kensington Maple

Solar Oak

Fonthill Pear

Wild Cherry

Windsor Mahogany

Shaker Cherry

Talas Cherry

Aria

Cocobala

Bed Specifications:

	Zenith 9200	Zenith 7200	Matrix 6200	Liberty
Adjustable Patient Surface Length:	80 - 84 inch	80 - 84 inch	80 - 84 inch	76 - 80 inch
Welded Sleep Deck Width:	35 inch	35 inch	35 or 42 inch	35 inch
Height Travel Range:	7 - 30 inch	7.95 - 30 inch	8.95 - 27 inch	9.95 - 24.4 inch
Safe Work Load: <i>includes resident, support surface, head/footboards, and accessories</i>	600 lb	550 lb	450 lb	450 lb
Maximum Patient Weight: <small>EVENLY DISTRIBUTED¹</small>	550 lb	500 lb	400 lb	400 lb
Actuators:	4 DC	4 DC	4 DC	3 DC

Options & Features:

Chair Position <small>(APS)</small> :	■	■	■	—
Advanced Positioning Pendant:	■	■	■	—
Pendant Attaches to Either Side of Bed:	■	■	■	—
Trendelenburg/Reverse Trendelenburg Positioning: <small>via APS staff control</small>	■	■	■	—
Slide-W-I-D-E™ Option Available: <small>expands to 35, 39, or 42 inches</small>	■	■	—	—
Patented Super-Efficient Transmission: <small>Improved speed and less strain on the motor</small>	■	—	—	—
Hi-Low Actuators are Unloaded in Low Position:	■	■	■	—
Embedded Staff Control: ²	■	■	■	—
Upgrade Hand Controls and/or Staff Controls Anytime: <small>interchangeable, mix and/or match standard and/or APS controls</small>	■	■	■	—
Splash Proof Electronics:	■	■	■	—
Intuitive Pedal Lock or 4 Locking Casters:	■	■	■	■ ²
Bed Locks and Rolls at Any Height: <small>reducing care-giver back strain</small>	■	■	■	■
Silent Operating, Sealed Ball Bearings: <small>Maintenance-Free</small>	■	■	■	—
Taller Casters: <small>approximately 7 inches of clearance</small>	—	—	—	■
Headboard and Footboard Mounting:	Tubular	Low Profile	Low Profile	Low Profile
Built-in Push Bar and Breakage Reduction:	■	—	—	—
Maximum Clearance: <small>overbed tables and lifts</small>	■	■	■	5.1 inch [†]
Can be Used as Transport Devices:	■	■	■	—
FDA Guideline Compliance: <small>for Bed Entrapment Risk Reduction³</small>	■	■	■	—
Silver Solutions Powder Coating: <small>process provides antimicrobial surface protection⁶</small>	■	■	■	—
Made in the USA by Skilled Craftsmen: <small>built to withstand trouble-free use for many years</small>	■	■	■	—

Accessories:

Edema Foot Ratchet: ² <small>elevates the foot end of the bed</small>	■	■	■	—
Plastic Pan Deck: ²	■	■	■	—
Backup Battery: ² <small>(external or mounted on board)</small>	■	■	■	—
Extension Kit: ² <small>(converts from 84 to 88 inch length)</small>	■	■	—	—
Pin-on Wide Kit: ²	■	■	■	—

Warranty:

Main Frame & Welds	15 years	15 years	15 years	10 years
Control Box & Actuator Motors	5 years	5 years	4 years	2 years
Hand Control Pendant & Cabling	3 years	3 years	3 years	1 year
Staff Control & Cabling	3 years	3 years	3 years	—
Head & Footboards	1 year	1 year	1 year	1 year

Product Page Link

--	--	--	--

■ Standard ■ Available ■ Optional — Not Applicable/Available

* options available at additional cost - accessories sold separately † fixed base ¹ maximum weight capacity - evenly distributed - refer to user manual ² anytime add-on or with purchase
³ boards constructed with 3DL material ⁴ Zenith and Matrix models ⁵ Hospital Bed System Dimensional and Assessment Guidance to Reduce Entrapment* details at www.fda.gov
⁶ inhibits the growth of bacterial odors, resists stains caused by mildew and/or mold, preserves the cured coating during the service life of the bed, and inhibits the growth of organisms that can detract from the appearance and coating functionality ⁷ Fragala G, Perry B, Fragala M, Examining bed width as a contributor to risk of falls from bed in long term care. *Annals of Long-Term Care: Clinical Care and Aging*, 2012; 20(6):35-38 ⁸ GF Health Products, Inc. makes no independent representations or empirical statements about Silver Solutions. This information is reprinted from: TCI Powder Coatings, EPA Guidance for EPA Registration 40 CFR 152.25
⁹ Wireframe retainers are standard on all beds except Slide-W-I-D-E® models which use plastic corner retainers