





*The* **LIBERTY BED**™

 **Basic**  
**American**™

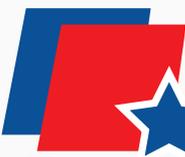

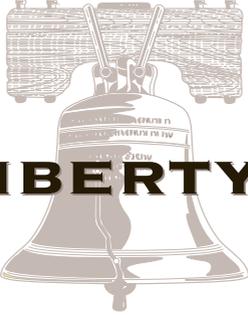

# The LIBERTY BED™

The Liberty Bed provides you with the key features you want at a very affordable price point.

Strong, no-sag grid deck sleep surface with sturdy cross braces.  
**Safe work load: 450 pounds\***

- **Height travel range: 9" to 24.4".** With a standard 6" mattress, this positions the bed at an ideal height for the caregiver and for transfers.



- **One bed can be used for 76" and 80" lengths** using the adjustable brackets to meet space limitations and improve inventory management. No extension kit is necessary.




- **Resident friendly, quiet, smooth operation** keeps bed level and balanced at all times using three DC motors.
- **360 degree mobility** at ANY height – can move from room to room at any height, and move over thresholds at even the lowest position.
- **Mattress retainer** is standard.
- **Reinforced, robust base** is designed for support yet provides sufficient underbed clearance, ideal for using overbed tables and patient lifts.
- **Easy assembly** requires only a screwdriver.



**Keyed Wall Saver is standard** – easily moves over thresholds during transport, yet automatically falls back into place to protect your walls. Wallsaver adjusts easily when modifying bed from 76" to 80" and when using a trapeze.

### Pendant Hand Control

Easy operation contoured Hand Control Pendant is shock-resistant, water-proof and durable. Holster conveniently secures to bedrails and assist devices for easy access.

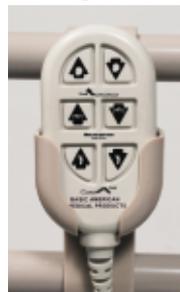



### Highest Safety and Quality Standards

Built to conform to UL 60601-1 and IEC601-2-38, CAN/CSA C22.2 No. 601.1 the Liberty Bed offers the following ensuring the utmost resident and staff safety:

- No pinch points within 200mm of bed perimeter.
- Rolling clearance of 20mm over thresholds.
- Splash-proof electronics.
- Load capacity and tip-stability testing.

### SilverSolutions within the powder-coating provides antimicrobial surface protection.

SilverSolutions inhibits the growth of organisms that can detract from the appearance and coating functionality:

- Inhibits the growth of bacterial odors
- Resists stains caused by mildew
- Resists stains caused by mold

SilverSolutions preserves the cured coating during the service life of the bed.

*GF Health Products, Inc. makes no independent representations or empirical statements about SilverSolutions. This information is reprinted from: TCI Powder Coatings, EPA Guidance for EPA Registration 40 CFR 152.25*

## ACCESSORIES & OPTIONS



**Optional Trapeze Adapter and Trapeze** attaches easily to headboard to position trapeze within easy reach of residents.  
ZA83000 (Adapter for use with ZA78100 Trapeze Bar)



**Optional Extension kit** Extends sleeping surface from 80" to 84".  
LBT83300

**Optional Liberty Lock** One-step red/green pedal locking system maximizes security. Locked position can be seen from a distance. The Liberty Brake can be added as a retrofit.



LBTL7174  
With Pedal Lock

**Optional Half Length Assist Device and Assist Bar** Half length assist device and assist bar provide a sturdy and secure handhold to assist residents in and out of bed. Assist devices are removable without tools and meet FDA Entrapment Guidelines.



ZA80300



ZA82800  
(head end, pair)  
ZA82900  
(foot end, pair)



ZA85500



**Four Locking Casters** secure the bed at any height.  
LBT7174 With 4 Locking Casters



**Basic American  
Liberty Bed™ Specifications**

Additional Features Include:  
Electrostatic powder-coated antimicrobial bed surface provides longer lasting, durable finish

- LBT7174** The Liberty Bed with four locking casters
- LBTL7174** The Liberty Bed with pedal lock

**Specifications:**

Height Travel Range .....	9.9" to 24.4"
Safe Work Load* .....	450 lbs.
Mattress Deck Length.....	76" or 80"
Mattress Deck Width.....	35"
Overall Bed Length.....	86" with wallsaver
Overall Bed Width.....	39.75" with rails/assist devices
Maximum back angle .....	70°
Maximum knee angle .....	30°

*Proudly Made in the U.S.A.* 

\*including resident, support surface, head/foot boards, and accessories.

For more information on the Liberty Bed, custom options and available accessories, please visit us online at **www.grahamfield.com** or call Basic American today.



**1.678.291.3200 Fax 1.770.368.4954**

Information contained herein is subject to change.  
The most current and complete product information can be found on our website  
**www.grahamfield.com**



© 2016, GF Health Products, Inc. All Rights Reserved.  
Basic American Medical Products is a trademark of GF Health Products, Inc.  
GF Health Products, Inc. is an ISO 13485 Certified Company