Basic Zenith 7100






- ★ Smart Technology
- ★ Digital Electronics enable 'plug and play' capability to add features at any time
 - USB Power Supply for cell phone charging
 - Underbed Light for resident safety
- ★ Interchangeable Electronics


Key Features:

- Energy-saving electronics consume up to 66% less power than electronics using a transformer
- Sealed ball bearings in the joints provide maintenance-free continuous, silent operation
- Intuitive 2-pedal locking system centrally located, easily visible from a distance to ensure the bed is locked, and is accessible at all heights and widths - even the lowest position
- US Patent awarded design of high/low mechanism inside of the bottom of the mainframe tubes reduces the chance of fluids getting inside and harboring bacteria
- Synchronized high-low motors keep the bed level for proper therapeutic positioning
- Industry-leading electronics ensure lower cost of ownership
- Tool-less assembly
- SilverSolutions within the powder-coating provides antimicrobial surface protection and deters the growth of stain and odor causing bacteria
- · Certified for use as a transport device



High Position Height Travel Range: 30"

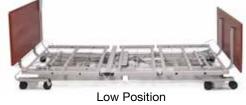

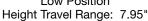



Chair Position



Trendelenburg Position



Reverse Trendelenburg

Advanced Positioning, including Cardiac Chair and Trendeleburg/Reverse Trendeleburg Safe Work Load: 500 lbs.

Specifications	
Overall Bed Length (with boards and wallsaver) 76" 80"	82" 86"
Sleep Deck Width	35"
Bed Width Including Boards	36"
Maximum Height (floor to top of mattress support deck)	30"
Minimum Height (floor to top of mattress support deck)	7.95"
Maximum Head / Back Deck Angle	68°
Maximum Knee / Foot Deck Angle	25°
Maximum Trendelenburg/Reverse Trendelenburg Angle (APS models)	15°
Maximum Safe Working Load WITH WEIGHT EVENLY DISTRIBUTED - includes bedding, resident / patient, support surface, and accessories	500 lb (226.8 kg)
Bed Mass (standard deck, without assist devices or boards)	204 lb (92.95kg)

Zenith710	00 with Standard Functionality	
ZZ79174	76/80"L x 35"W grid deck, intuitive 2-pedal locking mechanism	
ZG79174	76/80"L x 35"W grid deck, 4 locking and 4 non-locking casters	
Zenith7100 with Standard Functionality with SLIDE-W-I-D-E Deck		
ZZ79674	76/80"L x 35"W with expandable/retractable sleep deck, intuitive 2-pedal locking mechanism	
Zenith7100 with Advanced Positioning		
ZZA79174	76/80"L x 35"W grid deck with advanced positioning, intuitive 2-pedal locking mechanism	
ZAG79174	76/80"L x 35"W grid deck with advanced positioning, 4 locking and 4 non-locking casters	
Zenith7100 with Advance Positioning with SLIDE-W-I-D-E Deck		
ZZA79674	76/80"L x 35"W with advanced positioning, expandable/ retractable sleep deck, intuitive 2-pedal locking mechanism	

Optional Accessories





Counter-rotating Assist Device ZA84200 US Patent-awarded counter-rotating assist device





Pivot Assist Device za90300



Fixed Assist Device ZA85000





* Basic

American



Length Extension Kit **ZA90500** 4" extension ZA90400 8" extension

Warranty:

Main Frame and welds:	Fifteen years
Control Box and Actuator Motors:	Five years
Hand Control Pendant, Staff Control, and Cabling:	Three years
Headboard and Footboard:	One year
All other durable components not listed above:	Two years



Expandable/Retractable Sleep Deck

Patented industry-leading Slide-W-I-D-E deck option provides built-in adjustability to change the bed width from 35", 39", and 42"







Underbed Light clearly illuminates where the resident gets in and out of the bed ZA90010

USB Power Supply eliminates using a separate wall outlet for cell phone and tablet charging ZA90020



Battery Backup 999-0831-003 Charger 999-0711-002



Trapeze Adapter (mount) ŽA79000



Trapeze (bar chain, triangle) **ŽA78100**



Board Mounted IV Socket ZA19000



Embedded staff control with safety-inspired two-stage lockouts enables caregivers to lock-out individual features to the resident



Zenith/Matrix Bed Cart ZA89900



Standard 999-0900-301SP



Advanced Positioning 999-0900-305SP



Advanced Positioning (chair and T/RT) ZA831306

Information contained herein is subject to change. The most current and complete product information can be found on our website. www.grahamfield.com









