

## LIFT SLINGS

# SURELIFT™ UNIVERSAL SLINGS


## WITH FULL HEAD SUPPORT

Universal slings (also called divided leg slings) are the easiest sling to use while transferring a resident from a seated position. These can also be used for lifting from a supine position, making them suitable for most lifting situations. They are an easy-to-fit general purpose sling designed for the majority of users. They can be used for transfers from a bed to a chair, wheelchair, or toilet, chair to chair, or from the floor to a bed or chair.

Universal slings should be used for residents that are totally or partially dependent, non-weight bearing, or bariatric. Those who have limited head control should use a universal sling with head support. The large opening accommodates toileting. Mesh slings should be used for bathing. Lumex Universal slings have four or six connection points and can be used with 2-, 4-, or 6-point spreader bars.

### Features

- Provides excellent patient comfort and security
- Reinforced back for support
- Non-mesh slings have padding on legs and back
- Height-Adjustable straps
- Easy to use
- Machine washable
- Six point connectivity is designed to enhance positioning and safety
- Designed to meet the requirements of HCPCS Code: E0621
- Six (6) month limited warranty

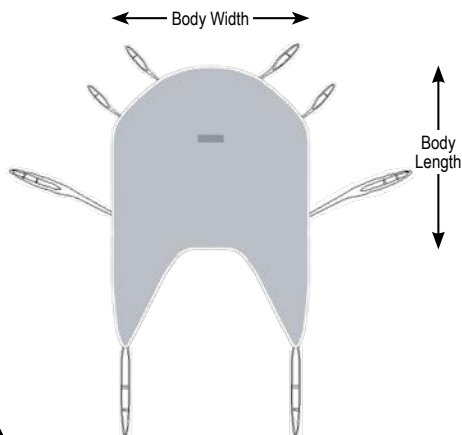

Size	Material	body Width		body Length		Weight Capacity EVENLY DISTRIBUTED		UOM
		inch	cm	inch	cm	lb	Kg	
UPH862	Medium	40	101.6	33.5	85.09	400	181.4	1 ea
UPH863	Large	43.5	110.49	40.5	102.87	400	181.4	1 ea
UPH864	XX-Large	46.75	118.75	44.25	112.4	600	272.2	1 ea
UMH832	Medium	40.5	102.87	34.5	87.63	400	181.4	1 ea
UMH833	Large	45	114.3	40.5	102.87	400	181.4	1 ea



GF Health Products, Inc. is an  
intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.  
The most current and complete product information can be found on our website.

[www.grahamfield.com](http://www.grahamfield.com)



© 2010, GF Health Products, Inc. All Rights Reserved.  
Lumex is a registered trademark of GF Health Products, Inc.

graham  
field

GF1000084RevC23